
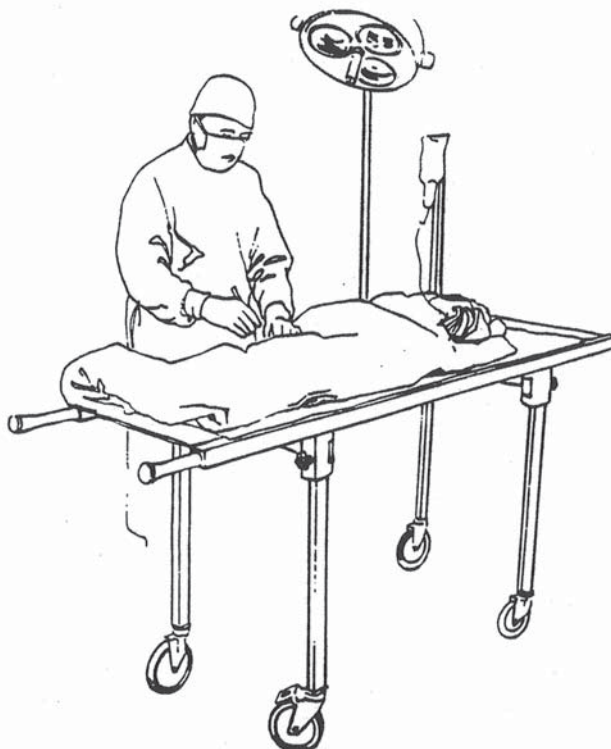
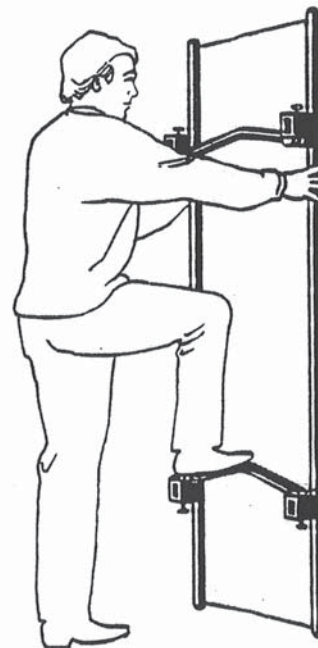


THE MULTI PURPOSE

# Norhosp Stretcher

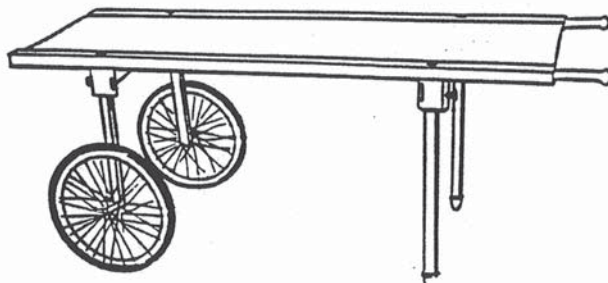

## The Multi Purpose NorHosp Stretcher

In all types of relief and emergency operations, it is of utmost importance that the available equipment is dependable, easy to use, and – above all – flexible.

With this goal in mind, **The Multi Purpose Norhosp Stretcher** was developed by Mediq Norge AS in cooperation with technologically advanced Norwegian companies and the Norwegian Defence.

This collapsible stretcher is made of aluminium sections and tarpaulin material which allows X-rays to pass through it, and the hand grips are moulded from extremely strong nylon.

**The Multi Purpose Norhosp Stretcher** can be equipped with stand poles – long, short or with wheels. The poles can be put through the side arms of the stretcher at whatever height desired. The stretcher can be placed perfectly level even rough ground, and several stretchers can be stacked on top of one another in order to accommodate maximum patients in a small area.

**The Multi Purpose Norhosp Stretcher** is the perfect solution for use in field dressing stations, at surgical hospitals and health clinics. The stretcher can be used for transportation, as a simple operating table, or as a hospital bed in a nursing unit. The stretcher is standard equipment in the NorHosp surgical field hospital's Nursing unit, and has frequently been used during operations worldwide for many years.

With the handgrips pushed back into the side arms, the total length of the standard **Multi Purpose Norhosp Stretcher** is 192cm, and will fit in all field vehicles, e.g. the Mercedes Gelandewagen, Land Rover etc. With the handgrips extended, the same stretcher fits into existing racks that are used in cars, boats, aircrafts and helicopters. Even in a small station wagon one can fit a stretcher by putting the front end on the seat of the front seat while supporting the rear end on its support stand poles.

### SPECIFICATIONS

Total length	220cm
Length of patient area	190cm
Width	57cm
Height	15cm
Stretcher weight	8,1kg
Each leg can tolerate	450 kg

### OPTIONAL EQUIPMENT

Support stand pole	60cm
Support stand pole	35cm
Support stand pole	w/wheel
Support stand pole	w/wheel and brake

Load tests have been conducted to check strength of material, welding methods etc. These are available upon request.